
Soprano Titanium

Stronger, Faster, Smarter

Hair Removal, **Reinvented.**



Introducing Soprano Titanium

A breakthrough hair removal platform that combines ultimate functionality with unparalleled comfort. Thanks to its large 4cm² spot size and advanced cooling system, treatments are now much faster and pain free – offering a significantly improved solution for patients and practitioners alike. Soprano Titanium leverages the benefits of three combined wavelengths, with an improved patient experience and a business-oriented approach, creating a unique and result-driven new solution in the world of professional hair removal.



The Titanium's System Highlights

✓ Exceptionally Large Spot Size

40% faster!

Soprano Titanium features TRIO 4cm², a unique and efficient new applicator with an especially large spot size of 4cm², capable of working on a grid of 600cm square. This breakthrough applicator makes treatments much faster and virtually painless, significantly cutting treatment times and ensuring a more comfortable experience for patients.

✓ Dual Connector

More Flexibility. Less Fuss.

Soprano Titanium is the first version of Soprano hair removal systems that includes two connectors, which enables two different applicators to be connected to the platform at all times. This feature allows for more flexibility, while eliminating time spent on switching applicators during treatments.

✓ ICEPlus⁺

Continuous Contact Cooling

ICE Plus is an advanced technology that continuously cools skin and facilitates range-controlled temperature throughout the entire treatment. The ICE Plus™ cold sapphire tip minimizes epidermal risks while maintaining heat within the dermis where hair follicles are treated. Temperature remains controlled and constant which contributes to the diode's efficiency. Patients can now enjoy a virtually painless, cool and comfortable experience throughout the session.

✓ 15" Android Screen

Friendly and Functional


Equipped with a new large sized 15" Android screen, with pre-set parameters, practitioners can enjoy a more intuitive, user friendly work experience.





Manage Your Clinic. Maximize Your Potential.

Smart Clinic is a pioneering, cloud-based businesses development tool that offers continuous access to live data through an application and online platform. It empowers better management processes simply and seamlessly, reflecting a real time status of platforms, applicators, and productivity, letting you access data and make adjustments without losing time. Smart Clinic drives optimization and standardization, empowering decisions based on genuine data and statistics – not guesswork. It also boosts marketing endeavors by enabling notifications of launches, new products, event invitations and more.

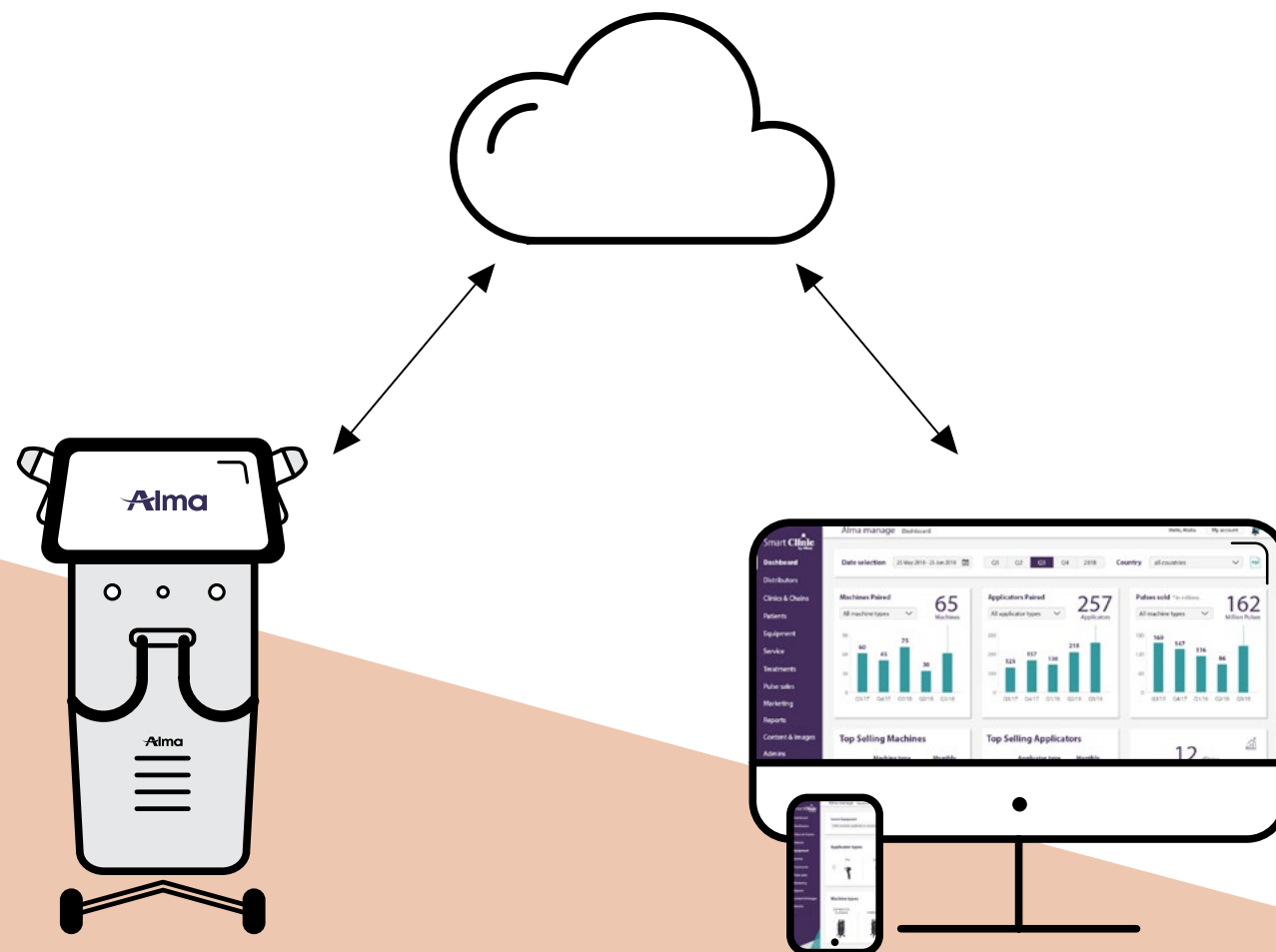

Never Miss a Beat

Smart Clinic is a multi-faceted solution that provides a complete picture of your business; whether it is a single clinic or a chain.

Treatment Management | Automatically document and log treatments, including areas treated, parameters used, and date of visit.

Platforms and Applicators | Lets you view the number of pulses, time operated, areas treated, number of platforms and applicators operated by each practitioner at any given time. This data is crucial for ongoing modifications of your clinic's equipment and for collecting operational statistics (per platforms, applicator and clinic), tracking practitioner performance and improving human resource management.

Maintenance and Service | Smart Clinic lets you open up a service call at any time through the application, sending you a troubleshoot guide so that problems can be solved when they arise. It offers the ability to track treatments and view statistics for each machine to prevent overuse and underuse, follow up on errors and receive ongoing maintenance alerts (water/filters). With Smart Clinic, you can prevent downtime, manage your resources more efficiently and boost productivity.





Synergetic Technologies. Superior Results.


SHR™ technique

SHR™ - Alma's proprietary hair removal innovation is the only clinically proven method of laser hair removal that is virtually painless, and it is the safest laser hair removal method for all skin types, including dark and tanned skin. Its unique gradual heating and high frequency pulse delivery method effectively damages the hair follicle and hinders re-growth, while preventing injury to the surrounding tissue. SHR™ protects the melanin of the skin using a unique method, allowing the procedure to be performed on all skin types anytime of the year as well as on tanned skin.

3D Technology



Featuring its exclusive 3D technology, Soprano Titanium combines the three most effective laser wavelengths into a single applicator, simultaneously targeting different tissue depths and anatomical structures within the hair follicle. By combining the absorption and penetration levels of different wavelengths, along with extended treatment coverage, increased comfort and low maintenance requirements, Soprano Titanium is the safest and most efficient hair removal treatment available today. 3D technology for one superior hair removal solution.

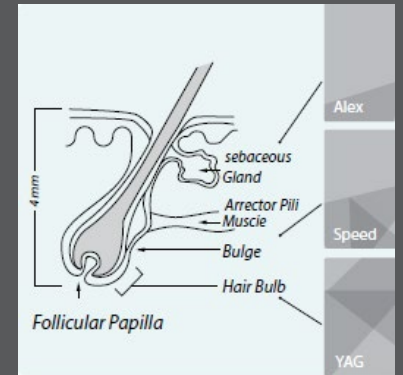

755nm

The 755nm wavelength, based on the Alex laser, offers a strong absorption by the chromophore melanin. This makes the 755nm an ideal wavelength in the removal of fine and light colored hair.

810nm

The 810nm wavelength is the gold standard for laser hair removal. It balances safety with efficacy by possessing good absorption in melanin and relatively deep penetration.

1064nm

The 1064nm wavelength is characterized by lower melanin absorption, making it a perfect solution for darker skin types. The 1064nm also offers deep penetration targeting the bulb and papilla of the hair follicle. With water and hemoglobin absorption, the 1064nm wavelength increases the thermal profile of the overall treatment, leading to highly effective hair removal, especially for deep hair follicles.

The Titanium's Applicators



TRIO 4cm²

3D and 40% faster!

This exclusive applicator is a valuable asset for any clinic. With its large 4cm² spot size, and integration of three wavelengths, it facilitates increased treatment coverage to significantly reduce treatment time and make the process easier and safer for patients. With TRIO 4cm², hair removal is quick, efficient and convenient.



TRIO - Three In One

This applicator enables the delivery of three highly effective wavelengths in one, making it suitable for all skin types. Covering a spot size of 2cm² it simultaneously targets different tissue depths and anatomical structures within the hair follicle, increasing treatment efficiency.



SPEED™ 810

Incorporating a highly effective and safe wavelength in a comfortable 2cm spot size.



ALEX 755

Powerful absorption, combined with continuous cooling, enables the treatment of light and fine hair for skin types I-V.



COMPACT 810

Incorporating a highly effective and safe wavelength, the 1cm spot size allows treating small areas comfortably.



FACIAL TIP

Easily treat areas that are normally hard to reach, including the ears, nostrils and glabella.

Platform Benefits



LARGE 4 cm² SPOT SIZE FOR QUICK TREATMENTS



ICE PLUS™ CONTINUOUS COOLING SYSTEM



TWO CONNECTORS



SMART CLINIC DATA CENTER



CLINICALLY PROVEN 3D TECHNOLOGY



INTUITIVE, USER FRIENDLY LARGE SCREEN (15" ANDROID) WITH PRE-SET PARAMETERS



I-VI SKIN TONES AND TANNED SKIN



NO CONSUMABLES

Alma Air

Understanding the nature of your practice, and the need for maximum ease and mobility, we developed Alma Air, an adjustable belt designed to carry Alma's applicators while you work. Alma Air helps lift the load off your back and make your work easy and effortless.

Its functional design eliminates the need to hold the umbilical cord while you work, so that you don't have to carry any excess weight. Convenient, comfortable, and made from a breathable material to prevent perspiration, Alma Air supports your back – as well as your clinic's efficiency.



“

I am very pleased with the results I am able to achieve with the Soprano Titanium. Besides for the high efficiency and safety record, which are the base for any hair removal treatment, there are additional features, such as Smart Clinic and the unique TRIO 4cm² applicator, which set it above and beyond any other system on the market today.

Dr. Pablo Naranjo MD, PhD, Medical Director of
Elite Laser Clinic & Laser Unit at NISA Hospital, Madrid, Spain

”

Soprano Titanium

Hair Removal, **Reinvented.**

Alma
For You. For Life.

Alma Lasers GmbH

Nordostpark 100-102
90411 Nuremberg, Germany
Tel. + 49 911 / 89 11 29-0
Fax + 49 911 / 89 11 29-99

Email: info@almalasers.com

www.almalasers.com

© 2022 Alma. All rights reserved.
Alma, its logo, Soprano Titanium, TRIO, ICE Plus, Smart Clinic,
SHR & Speed are trademarks and/or registered trademarks
of Alma in the United States and/or other countries.

PBSP28111801_05

CONNECT WITH ALMA

